



Superspeed centrifuges

Rotor guide for Sorvall LYNX Superspeed Centrifuge Series

Find the ideal fit for your application

With a comprehensive rotor portfolio

Thermo Scientific™ Rotor	Max. speed (rpm)	Max. RCF (x g) ¹	K Factor ²	Capacity: No. tubes x mL
Fixed angle				
Thermo Scientific™ Fiberlite™ F9-6x1000 LEX Rotor	9,000	17,568	2,886	6 x 1000
Thermo Scientific™ Fiberlite™ F10-4x1000 LEX Rotor	10,500	20,584	2,767	4 x 1000
Thermo Scientific™ Fiberlite™ F12-6x500 LEX Rotor	12,000	24,471	1,388	6 x 500
Thermo Scientific™ Fiberlite™ F14-6x250y Rotor	14,000	30,240	1,699	6 x 250
Thermo Scientific™ Fiberlite™ F14-14x50cy Rotor	14,000	33,746	798	14 x 50
Thermo Scientific™ Fiberlite™ F20-12x50 LEX Rotor	20,000	51,428	455	12 x 50
Thermo Scientific™ Fiberlite™ F21-8x50y Rotor	20,000	47,850	744	8 x 50
Thermo Scientific™ Fiberlite™ F23-48x1.5 Rotor	23,000	57,368	199	48 x 1.5
Thermo Scientific™ A21-24x15c Rotor	21,500	63,049 / 63,049 ³	379 / 558 ³	24 x 15
Thermo Scientific™ A22-24x16 Rotor	22,000	60,063 / 60,063 ³	377 / 547 ³	24 x 16
Thermo Scientific™ A23-6x100 Rotor	23,500	62,976	503	6 x 100
Thermo Scientific™ A27-6x50 Rotor	27,000	79,057	485	6 x 50
Thermo Scientific™ A27-8x50 Rotor	27,000	87,207	408	8 x 50
Thermo Scientific™ T29-8x50 Rotor	29,000	100,605	354	8 x 50
Swinging bucket				
Thermo Scientific™ BIOFlex™ HC Rotor	5,500	7,068	5,522	4 x 1000
Thermo Scientific™ BIOFlex™ HS Rotor	7,000	10,025	4,889	4 x 400
Thermo Scientific™ TH13-6x50 Rotor	13,100	30,314	1,503	6 x 50
Continuous flow/ zonal				
Thermo Scientific™ TCF-20 Continuous Flow Rotor	20,000	42,931	620	1350
Thermo Scientific™ TCF-20 Zonal Rotor	20,000	42,931	620	1350

¹ RCF calculated for maximum rotor radius.


² K factor calculated assuming maximum volume in a sealed tube.

³ Outer row value/inner row value.



Thermo Scientific™ Sorvall™ LYNX 6000 and 4000 Centrifuges

Thermo Scientific™ Fiberlite™ F9-6x1000 LEX Rotor

Rotor specifications	
Type	Fixed angle
Material	Carbon fiber composite
Capacity	6 x 1000 mL
Maximum speed	9,000 rpm
Tube dimensions Ø x L	98 x 195 mm
Angle	20°
K Factor	2,886
Net weight	15.7 kg/34.5 lbs
Warranty	15 years
Certifications	



	RCF	Radius
Maximum	17,568 x g	194 mm
Minimum	5,886 x g	77mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	9,000 rpm	N/A
Max. RCF	17,568 x g	N/A

Fiberlite F9-6x1000 LEX rotor package includes:

Description	Cat. No.
Fiberlite F9-6x1000 LEX Rotor, including the below items:	096-061075
1000 mL Thermo Scientific™ Fiberlite™ PPCO High Performance Bottle with Nylon Cap and PP Plug (qty. 6)	010-1491
Rotor lid assembly	099-061075
O-Ring refresher kit ¹	021-061075

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
1000.0	1000.0	Thermo Scientific™ Fiberlite™ PPCO Bottle ⁴	010-1491	2	6	9,000	17,568	98 x 195	Closure	Included	2	Nylon Cap w/ PP Plug
1000.0	1000.0	Thermo Scientific™ Fiberlite PC Bottle ⁴	010-1492	2	6	9,000	17,568	98 x 195	Closure	Included	2	Nylon Cap w/ PP Plug
1000.0	1000.0	Thermo Scientific™ Sorvall™ PPCO High Perf. Bottle ^{2, 4}	010-1456	2	6	9,000	17,568	98 x 195	Closure	75003511 Included	2	High Performance AI Cap w/ PP Plug
1000.0	1000.0	Sorvall PPCO High Perf. Bottle ^{3, 4}	010-1459	2	6	9,000	17,568	98 x 195	Closure	75003511 Included	2	High Performance AI Cap w/ PP Plug

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

¹ Includes 5 O-Rings and 3 O-Ring Grease Packets.

² Replacement PPCO bottle without closure, set of 2, 010-1455.

³ Replacement PC bottle without closure, set of 2, 010-1458.

⁴ Follow all chemical compatibility limitations.

Refer to manufacturer's recommendations for product performance information.

Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
500.0	400.0	Fiberlite PPCO Bottle	010-1493	6	6	9,000	17,568	70 x 160	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter	010-0145	2	1 place/adapter
500.0	400.0	Fiberlite PC Bottle	010-1494	6	6	9,000	17,568	70 x 160	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter	010-0145	2	1 place/adapter
500.0	450.0	Thermo Scientific™ Nalgene™ PPCO Bottle	3141-0500	24	6	7,950	13,700	70 x 160	Closure	Included	24	PP Sealing
									Adapter	010-0145	2	1 place/adapter
500.0	450.0	Thermo Scientific™ Nalgene™ PC Bottle	3140-0500	24	6	7,950	13,700	70 x 160	Closure	Included	24	PP Sealing
									Adapter	010-0145	2	1 place/adapter
250.0	250.0	Fiberlite PPCO Bottle	010-1495	6	6	9,000	17,568	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter	010-0150	2	1 place/adapter
250.0	250.0	Fiberlite PC Bottle	010-1496	6	6	9,000	17,568	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter	010-0150	2	1 place/adapter
250.0	250.0	Nalgene PPCO Bottle	3141-0250	36	6	9,000	17,568	61 x 133	Closure	Included	36	PP Sealing
									Adapter	010-0150	2	1 place/adapter
250.0	250.0	Nalgene PC Bottle	3140-0250	36	6	9,000	17,568	61 x 135	Closure	Included	36	PP Sealing
									Adapter	010-0150	2	1 place/adapter
250.0	-	Corning™ Disposable Conical Tube	-	-	6	9,000	17,568	-	Adapter	010-1096	2	1 place/adapter
175.0	175.0	Thermo Scientific™ Nalgene™ PPCO Conical Wide-Mouth Tube	3143-0175	36	6	9,000	17,568	61 x 144	Closure	Included	36	PP Sealing
									Adapter	010-1132	2	1 place/adapter
175.0	175.0	Thermo Scientific™ Nalgene™ PC Conical Wide-Mouth Tube	3144-0175	36	6	9,000	17,568	61 x 144	Closure	Included	36	PP Sealing
									Adapter	010-1132	2	1 place/adapter
100	100	Pyrex™ Tube	-	-	12	-	-	43 x 139	Adapter	010-1425	2	2 places/adapter
80.0	76.0	PPCO Oak Ridge Tube	010-1280	6	18	9,000	17,568	38 x 109	Closure	Included	6	PP Sealing
									Adapter	010-1093	2	3 places/adapter
85.0	81.0	Thermo Scientific™ Nalgene™ PC Oak Ridge Tube	3118-0085	100	18	9,000	17,568	38 x 109	Closure	Included	100	PP Screw Top
									Adapter	010-1093	2	3 places/adapter
50.0	50.0	Thermo Scientific™ Nunc™ Disposable Conical Tube	339653	25	18	9,000	17,568		Closure	Included	25	PP Sealing
									Adapter	010-0180	2	5 places/adapter
50.0	-	Corning™ Disposable Conical Tube	-	-	30	-	-	-	Adapter	010-0180	2	5 places/adapter
50.0	-	Falcon™ Disposable Conical Tube	-	-	30	-	-	-	Adapter	010-0180	2	5 places/adapter
50.0	46.0	Nalgene PC Oak Ridge Tube	3138-0050	50	30	9,000	17,568	29 x 115	Closure	Included	50	PP Sealing
									Adapter	010-1091	2	6 places/adapter


Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0	46.0	Nalgene PPCO Oak Ridge Tube	3139-0050	50	36	9,000	17,568	29 x 115	Closure	Included	50	PP Sealing
									Adapter	010-1091	2	6 places/adapter
50.0	46.0	Nalgene PSF Oak Ridge Tube	3137-0050	50	36	9,000	17,568	29 x 115	Closure	Included	50	PP Sealing
									Adapter	010-1091	2	6 places/adapter
50.0	46.0	Nalgene FEP Oak Ridge Tube	3114-0050	10	36	9,000	17,568	29 x 108	Closure	Included	10	ETFE Screw Cap
									Adapter	010-1091	2	6 places/adapter
30.0	30.0	Nalgene PC Oak Ridge Tube	3138-0030	50	42	9,000	17,568	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1095	2	7 places/adapter
30.0	30.0	Nalgene PPCO Oak Ridge Tube	3139-0030	50	42	9,000	17,568	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1095	2	7 places/adapter
30.0	30.0	Nalgene FEP Oak Ridge Tube	3114-0030	10	42	9,000	17,568	26 x 102	Closure	Included	10	ETFE Screw Cap
									Adapter	010-1095	2	7 places/adapter
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	84	9,000	17,568	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-1087	2	14 places/adapter
16.0	16.0	Nalgene PPCO Oak Ridge Tube	3139-0016	50	84	9,000	17,568	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-1087	2	14 places/adapter
15.0	-	Corning™ Disposable Conical Tube	-	-	72	-	-	-	Adapter	010-1079	2	12 places/adapter
15.0	-	Falcon™ Disposable Conical Tube	-	-	72	-	-	-	Adapter	010-1079	2	12 places/adapter
10.0	8.0	Nalgene PPCO Oak Ridge Tube	3139-0010	50	108	9,000	17,568	16 x 100	Closure	Included	50	PP Sealing
									Adapter	010-1307	2	18 places/adapter
10.0	8.0	Nalgene PC Oak Ridge Tube	3138-0010	50	108	9,000	17,568	16 x 100	Closure	Included	50	PP Sealing
									Adapter	010-1307	2	18 places/adapter
10.0	-	BD Vacutainer™ Tube	-	-	108	-	-	16 x 100	Adapter	010-1415	2	18 places/adapter
6.0	6.0	BD Vacutainer Tube	-	-	132	-	-	13 x 100	Adapter	010-1416	2	22 places/adapter
2.0	2.0	Filtration Tube and 1.5 mL Conical Tube	-	-	72	-	-	13 x 45	Adapter	010-1417	2	12 places/adapter
1.8-2.7	-	BD Vacutainer Tube	-	-	180	-	-	16 x 100	Adapter	010-1419	2	30 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ Fiberlite™ F10-4x1000 LEX Rotor

Rotor specifications

Type	Fixed angle
Material	Carbon fiber composite
Capacity	4 x 1000 mL
Maximum speed	10,500 rpm
Tube dimensions Ø x L	98 x 195 mm
Angle	20°
K Factor	2,767
Net weight	10.9 kg/24 lbs
Warranty	15 years
Certifications	



	RCF	Radius
Maximum	20,584 x g	167 mm
Minimum	3,205 x g	50 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	10,500 rpm	10,500 rpm
Max. RCF	20,584 x g	20,584 x g

Fiberlite F10-4x1000 LEX rotor package includes:

Description	Cat. No.
Fiberlite F10-4x1000 LEX Rotor, including the below items:	096-041075
1000 mL Fiberlite PPCO high performance bottle with nylon cap and PP plug (qty. 4)	010-1491
Rotor lid assembly	099-041075
O-Ring refresher kit ¹	021-041075

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./ pkg.	Tubes/ rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./ pkg.	Description
1000.0	1000.0	Fiberlite PPCO Bottle ⁴	010-1491	2	4	10,500	20,584	98 x 195	Closure	Included	2	Nylon Cap w/ PP Plug
1000.0	1000.0	Fiberlite PC Bottle ⁴	010-1492	2	4	10,500	20,584	98 x 195	Closure	Included	2	Nylon Cap w/ PP Plug
1000.0	1000.0	Sorvall PPCO High Perf. Bottle ⁴	010-1456	2	4	10,500	20,584	98 x 195	Closure	75003511 Included	2	High Performance AI Cap w/ PP Plug
1000.0	1000.0	Sorvall PC High Perf. Bottle ^{3,4}	010-1459	2	4	10,500	20,584	98 x 195	Closure	75003511 Included	2	High Performance AI Cap w/ PP Plug
500.0	400.0	Fiberlite PPCO Bottle	010-1493	6	6	10,500	18,859	70 x 160	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter		2	1 place/adapter

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

¹ Includes 5 O-Rings and 3 O-Ring Grease Packets.

² Replacement PPCO bottle without closure, set of 2, 010-1455.

³ Replacement PC bottle without closure, set of 2, 010-1458.

⁴ Follow all chemical compatibility limitations.

Refer to manufacturer's recommendations for product performance information.

Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
500.0	400.0	Fiberlite PC Bottle	010-1494	6	6	10,500	18,859	70 x 160	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter	010-0145	2	1 place/adapter
500.0	450.0	Nalgene PPCO Bottle	3141-0500	24	4	8,550	13,700	70 x 160	AI Cap Closure	Included	24	PP Sealing
									Adapter	010-0145	2	1 place/adapter
500.0	450.0	Nalgene PC Bottle	3140-0500	24	4	8,550	13,700	70 x 160	Closure	Included	24	PP Sealing
									Adapter	010-0145	2	1 place/adapter
250.0	250.0	Fiberlite PPCO Bottle	010-1495	6	4	10,500	18,366	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter	010-0150	2	1 place/adapter
250.0	250.0	Fiberlite PC Bottle	010-1496	6	4	10,500	18,366	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter	010-0150	2	1 place/adapter
250.0	250.0	Nalgene PPCO Bottle	3141-0250	36	4	10,500	18,366	61 x 133	Closure	Included	36	PP Sealing
									Adapter	010-0150	2	1 place/adapter
250.0	250.0	Nalgene PC Bottle	3140-0250	36	4	10,500	18,366	61 x 133	Closure	Included	36	PP Sealing
									Adapter	010-0150	2	1 place/adapter
250.0	250.0	Corning™ Disposable Conical Tube	-	-	4	10,500	18,366	-	Adapter	010-1096	2	1 place/adapter
175.0	-	Nalgene PPCO Conical Wide-Mouth Bottle	3143-0175	36	4	10,500	15,777	61 x 144	Closure	Included	36	PP Sealing
									Adapter	010-1132	2	1 place/adapter
175.0	-	Nalgene PC Conical Wide-Mouth Tube	3144-0175	36	4	10,500	15,777	61 x 144	Closure	Included	36	PP Sealing
									Adapter	010-1132	2	1 place/adapter
100	100	Pyrex™ Tube	-	-	20	-	-	43 x 139	Adapter	010-1425	2	2 places/adapter
85.0	81.0	Nalgene PC Oak Ridge Tube	3118-0085	100	12	10,500	19,352	38 x 109	Closure	Included	100	PP Sealing
									Adapter	010-1093	2	3 places/adapter
80.0	76.0	PP Oak Ridge Tube	010-1280	6	12	10,500	19,352	38 x 109	Closure	Included	6	PP Sealing
									Adapter	010-1093	2	3 places/adapter
50.0	-	Nunc Disposable Conical Tube	339653	25	20	-	-	30 x 115	Closure	Included	25	PP Sealing
									Adapter	010-0180	2	5 places/adapter
50.0	-	Corning™ Disposable Conical Tube	-	-	20	-	-	-	Adapter	010-0180	2	5 places/adapter
50.0	-	Falcon™ Disposable Conical Tube	-	-	20	-	-	-	Adapter	010-0180	2	5 places/adapter
50.0	46.0	Nalgene PC Oak Ridge Tube	3138-0050	50	28	10,500	18,119	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-1091	2	7 places/adapter


Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./Pkg.	Description
50.0	46.0	Nalgene PPCO Oak Ridge Tube	3139-0050	50	28	10,500	18,119	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-1091	2	7 places/adapter
50.0	46.0	Nalgene PSF Oak Ridge Tube	3137-0050	50	28	10,500	18,119	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-1091	2	7 places/adapter
50.0	46.0	Nalgene FEP Oak Ridge Tube	3114-0050	10	28	10,500	18,119	29 x 108	Closure	Included	10	ETFE Screw Cap
									Adapter	010-1091	2	7 places/adapter
30.0	30.0	Nalgene PC Oak Ridge Tube	3138-0030	50	28	10,500	18,119	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1095	2	7 places/adapter
30.0	30.0	Nalgene PPCO Oak Ridge Tube	3139-0030	50	28	10,500	18,119	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1095	2	7 places/adapter
30.0	30.0	Nalgene FEP Oak Ridge Tube	3114-0030	10	28	10,500	18,119	26 x 102	Closure	Included	10	ETFE Screw Cap
									Adapter	010-1095	2	7 places/adapter
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	60	10,500	18,119	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-1087	2	15 places/adapter
16.0	16.0	Nalgene PPCO Oak Ridge Tube	3139-0016	50	60	10,500	18,119	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-1087	2	15 places/adapter
15.0	-	Corning™ Disposable Conical Tube	-	-	48	-	-	-	Adapter	010-1079	2	12 places/adapter
15.0	-	Falcon™ Disposable Conical Tube	-	-	48	-	-	-	Adapter	010-1079	2	12 places/adapter
10.0	8.0	Nalgene PPCO Oak Ridge Tube	3139-0010	50	72	10,500	20,584	16 x 100	Closure	Included	50	PP Sealing
									Adapter	010-1307	2	18 places/adapter
10.0	8.0	Nalgene PC Oak Ridge Tube	3138-0010	50	72	10,500	20,584	16 x 100	Closure	Included	50	PP Sealing
									Adapter	010-1307	2	18 places/adapter
10.0	-	BD Vacutainer™ Tube	-	-	72	-	-	16 x 100	Adapter	010-1415	2	18 places/adapter
6.0	6.0	BD Vacutainer Tube	-	-	88	-	-	13 x 100	Adapter	010-1416	2	22 places/adapter
2.0	2.0	Filtration Tube and 1.5 mL Conical Tube	-	-	48	-	-	13 x 45	Adapter	010-1417	2	12 places/adapter
1.8-2.7	-	BD Vacutainer Tube	-	-	120	-	-	10 x 64	Adapter	010-1419	2	30 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ Fiberlite™ F12-6x500 LEX Rotor

Rotor specifications

Type	Fixed angle
Material	Carbon fiber composite
Capacity	6 x 500 mL
Maximum speed	12,000 rpm
Tube dimensions Ø x L	70 x 160 mm
Angle	20°
K Factor	1,388
Net weight	9.1 kg/20 lbs
Warranty	15 years
Certifications	



	RCF	Radius
Maximum	24,471 x g	152 mm
Minimum	6,601 x g	69 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	12,000 rpm	12,000 rpm
Max. RCF	24,471 x g	24,471 x g

Fiberlite F12-6x500 LEX rotor package includes:

Description	Cat. No.
Fiberlite F12-6x500 LEX rotor, including the below items:	096-062375
500 mL Fiberlite PPCO bottle with PPGF cap and PP plug (qty. 6)	010-1493
Rotor lid assembly	099-062375
O-Ring refresher kit ¹	021-062375

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
500.0	400.0	Fiberlite PPCO Bottle	010-1493	6	6	12,000	24,471	70 x 160	Closure	Included	6	PPGF Cap w/ PP Plug
500.0	400.0	Fiberlite PC Bottle	010-1494	6	6	12,000	24,471	70 x 160	Closure	Included	6	PPGF Cap w/ PP Plug
500.0	450.0	Nalgene PPCO Bottle	3141-0500	24	6	9,000	13,700	70 x 160	Closure	Included	24	PP Sealing
500.0	450.0	Nalgene PC Bottle	3140-0500	24	6	9,000	13,700	70 x 160	Closure	Included	24	PP Sealing
250.0	250.0	Fiberlite PPCO Bottle	010-1495	6	6	12,000	22,260	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter	010-0151	2	1 place/adapter

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

¹ Includes 5 O-Rings and 3 O-Ring Grease Packets.

Refer to manufacturer's recommendations for product performance information.

Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
250.0	250.0	Fiberlite PC Bottle	010-1496	6	6	12,000	22,260	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
									Adapter	010-0151	2	1 place/adapter
250.0	250.0	Nalgene PPCO Bottle	3141-0250	36	6	12,000	22,260	61 x 124	Closure	Included	36	PP Sealing
									Adapter	010-0151	2	1 place/adapter
250.0	250.0	Nalgene PC Bottle	3140-0250	36	6	12,000	22,260	61 x 124	Closure	Included	36	PP Sealing
									Adapter	010-0151	2	1 place/adapter
250.0	-	Corning™ Disposable Conical Bottle	-	-	6	-	-	-	Adapter	010-1135	2	1 place/adapter
175.0	-	Nalgene PPCO Conical Wide-Mouth Bottle	3143-0175	36	6	-	-	61 x 144	Closure	Included	36	PP Sealing
									Adapter	010-0152	2	1 place/adapter
175.0	-	Nalgene PC Conical Wide-Mouth Bottle	3144-0175	36	6	-	-	61 x 144	Closure	Included	36	PP Sealing
									Adapter	010-0152	2	1 place/adapter
80.0	73.0	PC Oak Ridge Tube	010-0515	6	6	12,000	19,190	38 x 109	Closure	Included	6	PP Sealing
									Adapter	010-1114	2	1 place/adapter
80.0	76.0	PPCO Oak Ridge Tube	010-1280	6	6	12,000	19,190	38 x 109	Closure	Included	25	PP Screw Top
									Adapter	010-1114	2	1 place/adapter
50.0	-	Nunc Disposable Conical Tube	339653	25	6	-	-	30 x 115	Closure	Included	25	PP Sealing
									Adapter	010-1102	2	1 place/adapter
50.0	-	Corning™ Disposable Conical Tube	-	-	6	-	-	-	Adapter	010-1102	2	1 place/adapter
50.0	-	Falcon™ Disposable Conical Tube	-	-	6	-	-	-	Adapter	010-1102	2	1 place/adapter
50.0	46.0	Nalgene PC Oak Ridge Tube	3138-0050	50	12	12,000	22,580	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-1112	2	2 places/adapter
50.0	46.0	Nalgene PPCO Oak Ridge Tube	3139-0050	50	12	12,000	22,580	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-1112	2	2 places/adapter
50.0	46.0	Nalgene PSF Oak Ridge Tube	3137-0050	50	12	12,000	22,580	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-1112	2	2 places/adapter
50.0	46.0	Nalgene FEP Oak Ridge Tube	3114-0050	10	12	12,000	22,580	29 x 108	Closure	Included	10	ETFE Screw Cap
									Adapter	010-1112	2	2 places/adapter
30.0	30.0	Nalgene PC Oak Ridge Tube	3138-0030	50	18	12,000	20,640	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1115	2	3 places/adapter


Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
30.0	30.0	Nalgene PPCO Oak Ridge Tube	3139-0030	50	18	12,000	20,640	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1115	2	3 places/adapter
30.0	30.0	Nalgene FEP Oak Ridge Tube	3114-0030	10	18	12,000	20,640	26 x 102	Closure	Included	10	ETFE Screw Cap
									Adapter	010-1115	2	3 places/adapter
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	42	12,000	20,640	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-1105	2	7 places/adapter
16.0	16.0	Nalgene PPCO Oak Ridge Tube	3139-0016	50	42	12,000	20,640	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-1105	2	7 places/adapter
15.0	15.0	Corning™ Disposable Conical Tube	-	-	36	-	-	-	Adapter	010-1099	2	6 places/adapter
15.0	15.0	Falcon™ Disposable Conical Tube	-	-	36	-	-	-	Adapter	010-1099	2	6 places/adapter
10.0	8.0	Nalgene PPCO Oak Ridge Tube	3139-0010	50	42	12,000	20,640	16 x 100	Closure	Included	50	PP Sealing
									Adapter	010-1308	2	7 places/adapter
10.0	8.0	Nalgene PC Oak Ridge Tube	3138-0010	50	42	12,000	20,640	16 x 100	Closure	Included	50	PP Sealing
									Adapter	010-1308	2	7 places/adapter
10.0	-	BD Vacutainer™ Tube	-	-	42	-	-	16 x 100	Adapter	010-1103	2	7 places/adapter
3.0	-	BD Vacutainer Tube	-	-	84	-	-	10 x 64	Adapter	010-1137	2	14 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ Fiberlite™ F14-6x250y Rotor



Rotor specifications	
Type	Fixed angle
Material	Carbon fiber composite
Capacity	6 x 250 mL
Maximum speed	14,000 rpm
Tube dimensions Ø x L	62 x 135 mm
Angle	23°
K Factor	1,699
Net weight	8.0 kg/17.5 lbs
Warranty	15 years
Certifications	

	RCF	Radius
Maximum	30,240 x g	138 mm
Minimum	8,108 x g	37 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	14,000 rpm	14,000 rpm
Max. RCF	30,240 x g	30,240 x g

Fiberlite F14-6x250y rotor package includes:

Description	Cat. No.
Fiberlite F14-6x250y rotor, including the below items:	096-062075
250 mL Fiberlite PPCO bottle with PPGF cap and PP plug (qty. 6)	010-1495
Rotor lid assembly	099-062075
O-Ring refresher kit ¹	022-062075

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
250.0	250.0	Fiberlite PPCO Bottle	010-1495	6	6	14,000	30,240	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
250.0	250.0	Fiberlite PC Bottle	010-1496	6	6	14,000	30,240	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
250.0	250.0	Nalgene PPCO Bottle	3141-0250	36	6	13,350	27,500	61 x 124	Closure	Included	36	PP Sealing
250.0	250.0	Nalgene PC Bottle	3140-0250	36	6	13,350	27,500	61 x 124	Closure	Included	36	PP Sealing
85.0	81.0	Nalgene PC Oak Ridge Tube	3118-0085	100	6	14,000	30,240	38 x 106	Closure	Included	100	PP Screw Top
									Adapter	010-1119	2	1 place/adapter

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

¹ Includes 5 O-Rings and 3 O-Ring Grease Packets.

Refer to manufacturer's recommendations for product performance information.


Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
80.0	76.0	PPCO Oak Ridge Tube	010-1280	6	6	14,000	30,240	38 x 109	Closure	Included	6	PP Sealing
									Adapter	010-1119	2	1 place/adapter
50.0	50.0	Nunc Disposable Conical Tube	339653	25	6	9,000	12,497	30 x 115	Closure	Included	25	PP Sealing
									Adapter	75100136	2	1 place/adapter
50.0	-	Corning™ Disposable Conical Tube	-	-	6	9,000	12,497	-	Adapter	75100136	2	1 place/adapter
50.0	-	Falcon™ Disposable Conical Tube	-	-	6	9,000	12,497	-	Adapter	75100136	2	1 place/adapter
50.0	46.0	Nalgene PC Oak Ridge Screw Top Tube	3138-0050	50	6	14,000	30,240	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-0138	2	1 place/adapter
50.0	46.0	Nalgene PPCO Oak Ridge Tube	3139-0050	50	6	14,000	30,240	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-0138	2	1 place/adapter
50.0	46.0	Nalgene PSF Oak Ridge Tube	3137-0050	50	8	14,000	30,240	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-0138	2	1 place/adapter
50.0	46.0	Nalgene FEP Oak Ridge Screw Top Tube	3114-0050	10	6	14,000	30,240	29 x 108	Closure	Included	10	ETFE Screw Cap
									Adapter	010-0138	2	1 place/adapter
30.0	30.0	Nalgene PC Oak Ridge Screw Top Tube	3138-0030	50	12	14,000	30,240	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1072	2	2 places/adapter
30.0	30.0	Nalgene PPCO Oak Ridge Screw Top Tube	3139-0030	50	12	14,000	30,240	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1072	2	2 places/adapter
30.0	30.0	Nalgene FEP Oak Ridge Tube	3114-0030	10	12	14,000	30,240	26 x 102	Closure	Included	10	ETFE Screw Cap
									Adapter	010-1072	2	2 places/adapter
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	30	14,000	30,240	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-1074	2	5 places/adapter
16.0	16.0	Nalgene PPCO Oak Ridge Tube	3139-0016	50	30	14,000	30,240	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-1074	2	5 places/adapter
15.0	-	Corning™ Disposable Conical Tube	-	-	30	-	-	-	Adapter	75101073	2	5 places/adapter
15.0	15.0	Nunc Conical Tube	339650	25	30	-	-	-	Adapter	75101073	2	5 places/adapter
15.0	-	Falcon™ Disposable Conical Tube	-	-	30	-	-	-	Adapter	010-1410	2	5 places/adapter
10.0	8.0	Nalgene PPCO Oak Ridge Tube	3139-0010	50	42	14,000	30,240	16 x 100	Closure	Included	50	PP Sealing
									Adapter	010-1309	2	7 places/adapter
10.0	8.0	Nalgene PC Oak Ridge Tube	3138-0010	50	42	14,000	30,240	16 x 100	Closure	Included	50	PP Sealing
									Adapter	010-1309	2	7 places/adapter
10.0	-	BD Vacutainer™ Tube	-	-	42	-	-	16 x 100	Adapter	010-1117	2	7 places/adapter
3.0	-	BD Vacutainer Tube	-	-	60	-	-	10 x 64	Adapter	010-1138	2	7 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ Fiberlite™ F14-14x50cy Rotor

Rotor specifications

Type	Fixed angle
Material	Carbon fiber composite
Capacity	14 x 50 mL
Maximum speed	14,000 rpm
Tube dimensions Ø x L	29 x 115 mm
Angle	34°
K Factor	798
Net weight	7.7 kg/17 lbs
Warranty	15 years
Certifications	



	RCF	Radius
Maximum	33,746 x g	154 mm
Minimum	18,188 x g	83 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	14,000 rpm	13,000 rpm
Max. RCF	33,746 x g	29,097 x g

Fiberlite F14-14x50cy rotor package includes:

Description	Cat. No.
Fiberlite F14-14x50cy rotor, including the below items:	096-145075
Rotor lid assembly	099-145075
O-Ring refresher kit ¹	021-145075

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0	45.0	Corning™ Disposable Conical Tube ²	-	-	14	14,000	33,746	-	-	-	-	-
50.0	45.0	Falcon™ Disposable Conical Tube2	-	-	14	14,000	33,746	-	-	-	-	-
50.0	45.0	Nunc Disposable Conical Tube2	339653	25	14	8,690	13,000	30 x 115	Closure	Included	25	PP Sealing
50.0	-	Amicon™ Filtration Tube	-	-	14	5,000	4,304	-	-	-	-	-
50.0	46.0	Nalgene PC Oak Ridge Tube	3138-0050	50	4	14,000	33,746	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-0377	2	1 place/adaptor

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

¹ Includes 5 O-Rings and 3 O-Ring Grease Packets.

² Listed Max Speeds for Conical Disposable tubes may be higher than OEM tube specifications and may result in crazing and other cosmetic damage to the tube. Do not use these tubes more than once.

Refer to manufacturer's recommendations for product performance information.


Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0	46.0	Nalgene PPCO Oak Ridge Tube	3139-0050	50	14	14,000	33,527	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-0377	2	1 place/adapter
50.0	46.0	Nalgene PSF Oak Ridge Tube	3137-0050	50	14	14,000	33,527	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-0377	2	1 place/adapter
50.0	46.0	Nalgene FEP Oak Ridge Tube	3114-0050	10	14	14,000	33,527	29 x 108	Closure	Included	10	ETFE Screw Cap
									Adapter	010-0377	2	1 place/adapter
30.0	30.0	Nalgene PC Oak Ridge Tube	3138-0030	50	14	14,000	33,527	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1147	2	1 place/adapter
30.0	30.0	Nalgene PPCO Oak Ridge Tube	3139-0030	50	14	14,000	33,527	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-1147	2	1 place/adapter
30.0	30.0	Nalgene FEP Oak Ridge Tube	3114-0030	50	14	14,000	33,527	26 x 102	Closure	Included	50	ETFE Screw Cap
									Adapter	010-1147	2	1 place/adapter
16.0	16.0	Nalgene PC Oak Ridge Tube	3139-0016	50	14	14,000	33,527	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-0376	2	1 place/adapter
16.0	16.0	Nalgene PPCO Oak Ridge Tube	3138-0016	50	14	14,000	33,527	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-0376	2	1 place/adapter
15.0	-	Corning™ Disposable Conical Tube	-	-	14	14,000	33,527	-	Adapter	75100378	2	1 place/adapter
15.0	-	Falcon™ Disposable Conical Tube	-	-	14	-	-	-	Adapter	75100378	2	1 place/adapter
15.0	-	Amicon Filtration Devices	-	-	14	5,400	4,988	-	Adapter	010-1340	2	1 place/adapter
15.0	15.0	Nunc EZ Flip™ Conical Tube	362694	50	14	9,000	13,855	-	Adapter	010-1340	2	1 place/adapter
15.0	15.0	Nunc Conical Disposable Tube	339650	50	14	9,000	13,855	17 x 120	Closure	Included	25	PP Sealing
									Adapter	010-1340	2	1 place/adapter
10.0	8.0	Nalgene PPCO Oak Ridge Tube	3139-0010	50	14	14,000	33,527	16 x 100	Closure	Included	10	PP Sealing
									Adapter	010-1311 2	2	1 place/adapter
10.0	8.0	Nalgene PC Oak Ridge Tube	3138-0010	50	14	14,000	33,527	16 x 100	Closure	Included	10	PP Sealing
									Adapter	010-1311	2	1 place/adapter
10.0	-	BD Vacutainer™ Tube	-	-	14	-	-	16 x 100	Adapter	010-1124	2	1 place/adapter
5.0	-	5ml Eppendorf Microtube	-	-	14	12,000	25,000	-	Adapter	75005770	1	1 place/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ Fiberlite™ F20-12x50 LEX Rotor



Rotor specifications		
Type	Fixed angle	
Material	Carbon fiber composite	
Capacity	12 x 50 mL	
Maximum speed	20,000 rpm	
Tube dimensions Ø x L	29 x 115 mm	
Angle	25°	
K Factor	455	
Net weight	5.25 kg/11.5 lbs	
Warranty	15 years	
Certifications		
	RCF	Radius
Maximum	51,428 x g	115 mm
Minimum	25,043 x g	56 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	20,000 rpm	18,000 rpm
Max. RCF	51,428 x g	41,657 x g

Fiberlite F20-12x50 LEX rotor package includes:

Description	Cat. No.
Fiberlite F20-12x50 LEX rotor, including the below items:	096-124375
50 mL Nalgene PPCO Oak Ridge tube with sealing cap (qty. 12)	3139-0050
Rotor lid assembly	099-124375
O-Ring refresher kit ¹	021-124375

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0	46.0	Nalgene PC Oak Ridge Tube	3138-0050	50	12	20,000	51,428	29 x 115	Closure	Included	50	PP Sealing
50.0	46.0	Nalgene PPCO Oak Ridge Tube	3139-0050	50	12	20,000	51,428	29 x 115	Closure	Included	50	PP Sealing
50.0	46.0	Nalgene PSF Oak Ridge Tube	3137-0050	50	12	20,000	51,428	29 x 115	Closure	Included	50	PP Sealing
50.0	46.0	Nalgene FEP Oak Ridge Tube	3114-0050	10	12	20,000	51,428	29 x 108	Closure	Included	10	ETFE Screw Cap
30.0	30.0	Nalgene PC Oak Ridge Tube	3138-0030	50	12	20,000	51,428	29 x 108	Closure	Included	50	PP Sealing
									Adapter	010-0167	2	1 place/adaptor

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

¹ Includes 5 O-Rings and 3 O-Ring Grease Packets.

Refer to manufacturer's recommendations for product performance information.

Labware ordering information continued


Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
30.0	30.0	Nalgene PPCO Oak Ridge Tube	3139-0030	50	12	20,000	51,428	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-0167	2	1 place/adapter
30.0	30.0	Nalgene FEP Oak Ridge Tube	3114-0030	10	12	20,000	51,428	26 x 102	Closure	Included	10	ETFE Screw Cap
									Adapter	010-0167	2	1 place/adapter
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	12	20,000	51,428	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-0382	2	1 place/adapter
16.0	16.0	Nalgene PPCO Oak Ridge Tube	3139-0016	50	12	20,000	51,428	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-0382	2	1 place/adapter
15.0	-	Corning™ Disposable Conical Tube3	-	-	12	-	-	-	Closure	Included	-	PP Sealing
									Adapter	010-1123	2	1 place/adapter
15.0	-	Falcon™ Disposable Conical Tube3	-	-	12	-	-	-	Adapter	010-1123	2	1 place/adapter
14.0	11.0	PC Flanged Tube	03246	50	12	20,000	51,428	18 x 75	Closure	03269	50	PP Snap-on
									Adapter	75003025	2	1 place/adapter
10.0	8.0	Nalgene PPCO Oak Ridge Tube	3139-0010	50	12	20,000	51,428	16 x 82	Closure	Included	50	PP Sealing
									Adapter	010-1306	2	1 place/adapter
10.0	8.0	Nalgene PC Oak Ridge Tube	3138-0010	50	12	20,000	51,428	16 x 82	Closure	Included	50	PP Sealing
									Adapter	010-1306	2	1 place/adapter
10.0	-	BD Vacutainer™ Tube	-	-	12	-	-	16 x 100	Adapter	010-1068	2	1 place/adapter
4.0	3.0	PP Flanged Tube	03105	50	24	19,200	47,850	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 places/adapter
4.0	3.0	PC Flanged Tube	03104	50	24	19,200	47,850	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 places/adapter
3.0	-	BD Vacutainer Tube	-	-	12	-	-	10 x 64	Adapter	010-1128	2	1 place/adapter
1.5	1.5	Polyallomer Microtube	314352H01	100	36	16,600	35,776	11 x 40	Adapter	75003029	2	3 places/adapter
1.0	1.0	BD Microtainer™ Tube	-	-	36	-	-	8 x 48	Adapter	010-1127	2	3 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ Fiberlite™ F21-8x50y Rotor



Rotor specifications

Type	Fixed angle
Material	Carbon fiber composite
Capacity	8 x 50 mL
Maximum speed	20,000 rpm
Tube dimensions Ø x L	29 x 115 mm
Angle	34°
K Factor	744
Net weight	4.8 kg/10.5 lbs
Warranty	15 years
Certifications	

	RCF	Radius
Maximum	47,850 x g	107 mm
Minimum	14,758 x g	33 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	20,000 rpm	18,000 rpm
Max. RCF	47,850 x g	38,759 x g

Fiberlite F21-8x50y rotor package includes:

Description	Cat. No.
Fiberlite F21-8x50y rotor, including the below items:	096-084275
50 mL Nalgene PPCO Oak Ridge tube with sealing cap (qty. 8)	3139-0050
Rotor lid assembly	099-084275
O-Ring refresher kit ¹	021-084275

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0	46.0	Nalgene PC Oak Ridge Tube	3138-0050	50	8	20,000	47,850	29 x 115	Closure	Included	50	PP Sealing
50.0	46.0	Nalgene PPCO Oak Ridge Tube	3139-0050	50	8	20,000	47,850	29 x 115	Closure	Included	50	PP Sealing
50.0	46.0	Nalgene PSF Oak Ridge Tube	3137-0050	50	8	20,000	47,850	29 x 115	Closure	Included	50	PP Sealing
50.0	46.0	Nalgene FEP Oak Ridge Tube	3114-0050	10	8	20,000	47,850	29 x 108	Closure	Included	10	ETFE Screw Cap
30.0	30.0	Nalgene PC Oak Ridge Tube	3138-0030	50	8	20,000	47,850	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-0167	2	1 place/adapter

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

¹ Includes 5 O-Rings and 3 O-Ring Grease Packets.

Refer to manufacturer's recommendations for product performance information.

Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
30.0	30.0	Nalgene PPCO Oak Ridge Tube	3139-0030	50	8	20,000	47,850	26 x 102	Closure	Included	50	PP Sealing
									Adapter	010-0167	2	1 place/adapter
30.0	30.0	Nalgene FEP Oak Ridge Tube	3114-0030	10	8	20,000	47,850	26 x 102	Closure	Included	10	ETFE Screw Cap
									Adapter	010-0167	2	1 place/adapter
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	8	20,000	47,850	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-0382	2	1 place/adapter
16.0	16.0	Nalgene PPCO Oak Ridge Tube	3139-0016	50	8	20,000	47,850	18 x 107	Closure	Included	50	PP Sealing
									Adapter	010-0382	2	1 place/adapter
15.0	-	Corning™ Disposable Conical Tube ³	-	-	8	-	-	-	Adapter	010-1123	2	1 place/adapter
15.0	-	Falcon™ Disposable Conical Tube ³	-	-	8	-	-	-	Adapter	010-1123	2	1 place/adapter
14.0	11.0	PC Flanged Tube	03246	50	8	20,000	47,850	18 x 75	Closure	03269	50	PP Snap-on
									Adapter	75003025	2	1 place/adapter
10.0	8.0	Nalgene PPCO Oak Ridge Tube	3139-0010	50	8	20,000	47,850	16 x 82	Closure	Included	50	PP Sealing
									Adapter	010-1306	2	1 place/adapter
10.0	8.0	Nalgene PC Oak Ridge Tube	3138-0010	50	8	20,000	47,850	16 x 82	Closure	Included	50	PP Sealing
									Adapter	010-1306	2	1 place/adapter
10.0	-	BD Vacutainer™ Tube	-	-	8	-	-	16 x 100	Adapter	010-1068	2	1 place/adapter
4.0	3.0	PP Flanged Tube	03105	50	16	20,000	47,850	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 places/adapter
4.0	3.0	PC Flanged Tube	03104	50	16	20,000	47,850	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 places/adapter
3.0	-	BD Vacutainer Tube	-	-	8	-	-	10 x 64	Adapter	010-1128	2	1 place/adapter
1.0	-	BD Microtainer™ Tube	-	-	24	-	-	8 x 48	Adapter	010-1127	2	3 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ Fiberlite™ F23-48x1.5 Rotor

Rotor specifications

Type	Fixed angle
Material	Carbon fiber composite
Capacity	48 x 1.5 mL
Maximum speed	23,000 rpm
Tube dimensions Ø x L	11 x 40 mm
Angle	45°
K Factor	199
Net weight	4.54 kg/10.0 lbs
Warranty	15 years



	RCF	Radius
Maximum	57,368 x g	97 mm
Minimum	37,851 x g	64 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	23,000 rpm	18,500 rpm
Max. RCF	57,368 x g	37,116 x g

Fiberlite F23-48x1.5 rotor package includes:

Description	Cat. No.
Fiberlite F23-48x1.5 rotor, including the below items:	096-484075
1.5 mL high performance microtubes, pack of 100	314352H01
Rotor lid assembly	099-484075
O-Ring refresher kit ¹	021-484075

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
1.5	1.5	Polyallomer Microtube	314352H01	100	48	23,000	57,368	11 x 40	Closure	Included	-	Snap-on Seal

Refer to manufacturer's recommendations for product performance information.


¹ Includes 5 O-Rings and 3 O-Ring Grease Packets.

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ A21-24x15c Rotor



Rotor specifications

Type	Fixed angle
Material	Aluminum
Capacity	24 x 15 mL
Maximum speed	21,500 rpm
Tube dimensions Ø x L	121.5 x 16.5 mm
Angle	42/30°
K Factor	379
Net weight	8.7 kg /19.1 lbs
Warranty	10 years
Certifications	

	RCF	Radius
Maximum	63,049 x g (outer) 63,049 x g (inner)	122 mm (outer) 22 mm (inner)
Minimum	31,525 x g (outer) 22,739 x g (inner)	61 mm (outer) 44 mm (inner)

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	21,500 rpm	19,500 rpm
Max. RCF	63,049 x g	51,865 x g

A21-24x15c rotor package includes:

Description	Cat. No.
A21-24x15c rotor, including the below items:	75003004
Vacuum grease	76003500
Lubricant	75003786
Instruction manual	50136234
Retaining ring pliers	65614

Accessories:

Description	Cat. No.
Rotor stand	75003711
Replacement O-Ring	75007003

Labware ordering information continued


Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
15.0	15.0	Nunc PP Conical Disposable Tube ¹	339651	500	24	11,600	18,500	-	Closure	Included	-	-
15.0	15.0	Greiner™ PP Conical Disposable Tube ¹	-	-	24	21,500	63,049	-	Closure	Included	-	-
15.0	15.0	Greiner™ Conical Disposable Tube ¹	-	-	24	21,500	63,049	-	Closure	Included	-	-
15.0	15.0	Falcon™ PP Conical Disposable Tube ¹	-	-	24	21,500	63,049	-	Closure	Included	-	-
15.0	15.0	Corning™ PP Conical Disposable Tube ¹	-	-	24	21,500	63,049	-	Closure	Included	-	-
15.0	15.0	Corning™ PET Conical Disposable Tube ¹	-	-	24	21,500	63,049	-	Closure	Included	-	-
15.0	15.0	Sarstedt™ PP Conical Disposable Tube ¹	-	-	24	21,500	63,049	-	Closure	Included	-	-

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ A22-24x16 Rotor



Rotor specifications

Type	Fixed angle
Material	Aluminum
Capacity	24 x 16 mL
Maximum speed	22,000 rpm
Tube dimensions Ø x L	18 x 102 mm
Angle	28/42°
K Factor	377 (outer); 547 (inner)
Net weight	7.7 kg /17.0 lbs
Warranty	10 years
Certifications	

	RCF	Radius
Maximum	60,063 x g (outer) 60,063 x g (inner)	111 mm (outer) 111 mm (inner)
Minimum	29,220 x g (outer) 21,103 x g (inner)	54 mm (outer) 39 mm (inner)

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	22,000 rpm	22,000 rpm
Max. RCF	60,063 x g	60,063 x g

A22-24x16 rotor package includes:

Description	Cat. No.
A22-24x16 rotor, including the below items:	75003005
16 mL Nalgene PPCO Oak Ridge tube with sealing cap (qty. 24)	3139-0016
Vacuum grease	76003500
Lubricant for metal parts	75003786
Retaining ring pliers	65614
Instruction manual CD	50136234

Accessories:

Description	Cat. No.
Rotor Stand	75003711

Labware ordering information


Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
16.0	13.0	PP Flanged Tube	03244	50	24	22,000	60,063	18 x 100	Closure	03299	50	PP Snap-on
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	24	22,000	60,063	18 x 107	Closure	Included	-	PP Sealing
16.0	16.0	Nalgene PPCO Oak Ridge Tube	3139-0016	50	24	22,000	60,063	18 x 107	Closure	Included	-	PP Sealing
7.0	6.0	PC Flanged Tube	03120	50	24	22,000	60,063	13 x 100	Closure	03265	50	PP Snap-on
									Adapter	00416	1	1 place/adapter
7.0	6.0	PP Flanged Tube	03121	50	24	22,000	60,063	13 x 100	Closure	03265	50	PP Snap-on
									Adapter	00416	1	1 place/adapter
5.0	5.0	Pyrex™ Tube	03102	50	24	5,000	3,102	12 x 75	Adapter	00376	1	1 place/adapter
1.5	1.5	Polyallomer Microtube	314352H01	100	24	22,000	60,063	11 x 40	Adapter	00376	1	1 place/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ A23-6x100 Rotor



Rotor specifications

Type	Fixed angle
Material	Aluminum
Capacity	6 x 100 mL
Maximum speed	23,500 rpm
Tube dimensions Ø x L	38.0 x 105.0 mm
Angle	25°
K Factor	503
Net weight	7.2 kg /15.8 lbs
Warranty	10 years
Certifications	

	RCF	Radius
Maximum	62,976 x g	102 mm
Minimum	20,992 x g	34 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	23,500 rpm	21,100 rpm
Max. RCF	62,976 x g	50,290 x g

A22-24x16 rotor package includes:

Description	Cat. No.
A23-6x100 rotor, including the below items:	75003006
85 mL Nalgene PPCO Oak Ridge tube with sealing cap (qty. 6)	3118-0085
Vacuum grease	76003500
Lubricant for metal parts	75003786
Instruction manual CD	50136234
Retaining ring pliers	65614

Accessories:

Description	Cat. No.
Rotor Stand	75003711

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
94.0	70.0	PP Thick-Walled Tube	76002872	1	8	23,500	62,976	38 x 105	Closure	75001568	1	Sealing Cap
94.0	90.0	PC Screw Cap Tube	75002829	1	8	23,500	62,976	38 x 105	Closure	75001540	1	Sealing Cap
94.0	70.0	PC Screw Cap Tube	75002810	1	8	23,500	62,976	38 x 105	Closure	75001568	1	Sealing Cap

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

Refer to manufacturer's recommendations for product performance information.

Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
85.0	71.0	Nalgene PC Oak Ridge Tube	3118-0085	100	8	23,500	62,976	29 x 114	Closure	Included	100	PP Sealing Cap
50.0	37.0	PP Flanged Tube	03147	25	8	23,500	62,976	29 x 102	Closure	03268	25	PP Snap-on
									Adapter	75003102	1	1 place/adapter
50.0	36.0	PC Flanged Tube	03146	25	8	23,500	62,976	29 x 102	Closure	03268	25	PP Snap-on
									Adapter	75003102	1	1 place/adapter
50.0	46.0	Nalgene FEP Oak Ridge Tube	3114-0050	10	8	23,500	62,976	29 x 108	Closure	Included	10	FEP Sealing
									Adapter	75003102	1	1 place/adapter
50.0	43.0	Nalgene PC Oak Ridge Tube	3138-0050	50	8	23,500	62,976	29 x 107	Closure	Included	50	PP Sealing
									Adapter	75003102	1	1 place/adapter
50.0	43.0	Nalgene PP Oak Ridge Tube	3139-0050	50	8	23,500	62,976	29 x 107	Closure	Included	50	PP Sealing
									Adapter	75003102	1	1 place/adapter
50.0	45.0	Conical Tissue Culture Tube	-	-	6	17,000	32,956	30 x 121	Adapter	75003103	-	1 place/adapter
50.0	43.0	Round Bottom Tube	-	-	6	23,500	62,976	30 x 117	Adapter	75003102	-	1 place/adapter
38.0	35.0	Round Bottom Tube	-	-	6	23,500	62,976	26 x 110	Adapter	75003094	-	1 place/adapter
30.0	28.0	PC Oak Ridge Tube	314348	25	8	23,500	62,976	25 x 94	Closure	314347	1	Aluminum Sealing Assembly
									Adapter	75003027	1	1 place/adapter
									Adapter	75003102	-	1 place/adapter
									Tool	314353	1	Tube Extractor Tool
30.0	28.0	PP Oak Ridge Tube	314349	25	8	23,500	62,976	25 x 94	Closure	314347	1	Aluminum Sealing Assembly
									Adapter	75003027	1	1 place/adapter
									Adapter	75003102	1	1 place/adapter
									Tool	314353	1	Tube Extractor Tool
30.0	28.0	PP Tube	75002951	10	8	23,500	62,976	18 x 100	Closure	75002957	10	PP Stopper
									Adapter	75003102	1	1 place/adapter
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	6	23,500	62,976	18 x 123	Adapter	76002906	-	1 place/adapter
16.0	16.0	Nalgene PP Oak Ridge Tube	3139-0016	50	6	23,500	62,976	18 x 123	Adapter	76002906	-	1 place/adapter
15.0	15.0	Conical Tissue Culture Tube	-	-	6	17,000	32,956	17 x 123	Adapter	75003095	-	1 place/adapter
12.0	11.0	PC Flanged Tube	03115	50	8	23,500	62,976	16 x 100	Closure	03266	50	PP Non-sealing
									Adapter	00402	1	1 place/adapter
									Adapter	75003102	1	1 place/adapter
12.0	12.0	Round Bottom Tube	-	-	12	23,500	62,976	16 x 95	Adapter	75003093	-	2 places/adapter
10.0	9.1	PP Oak Ridge Tube	03929	50	8	23,500	62,976	16 x 80	Closure	03279	25	PP Sealing
									Adapter	75003024	2	1 place/adapter
									Adapter	75003102	1	1 place/adapter


Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
10.0	9.1	PC Oak Ridge Tube	03020	25	8	23,500	62,976	16 x 80	Closure	03279	25	PP Sealing
									Adapter	75003024	2	1 place/adapter
									Adapter	75003102	1	1 place/adapter
7.0	6.0	PC Flanged Tube	03120	50	8	23,500	62,976	13 x 100	Closure	03265	50	PP Snap-on
									Adapter	00473	1	1 place/adapter
									Adapter	75003102	1	1 place/adapter
6.5	6.5	Round Bottom Tube	-	-	12	23,500	62,976	13 x 114	Adapter	75003092	-	2 places/adapter
4.0	3.0	PC Flanged Tube	03104	50	16	23,500	62,976	11 x 75	Closure	03264	50	PP Snap-on
									Adapter	75003023	2	2 places/adapter
									Adapter	75003102	1	1 place/adapter
3.5	3.5	Round Bottom Tube	-	-	24	23,500	62,976	11 x 100	Adapter	75003091	-	4 places/adapter
1.5/2.0	-	Conical Microtube	-	-	24	23,500	62,976	11 x 40	Adapter	76002905	-	4 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ A27-6x50 Rotor

Rotor specifications

Type	Fixed angle
Material	Aluminum
Capacity	6 x 50 mL
Maximum speed	27,000 rpm
Tube dimensions Ø x L	29 x 108 mm
Angle	34°
K Factor	485
Net weight	5.58 kg/12.3 lbs
Warranty	10 years
Certifications	



	RCF	Radius
Maximum	79,057 x g	97 mm
Minimum	19,561 x g	24 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	27,000 rpm	24,000 rpm
Max. RCF	79,057 x g	62,465 x g

A27-6x50 rotor package includes:

Description	Cat. No.
A27-6x50 rotor, including the below items:	75003007
50 mL Nalgene PPCO Oak Ridge tube with sealing cap (qty. 6)	3139-0050
O-Ring replacement kit	75007006
Retaining ring pliers	65614
Instruction manual CD	50136234
Lubricant for metal parts	75003786
Vacuum grease	76003500

Accessories:

Description	Cat. No.
Ultracrimp sealing tool and crimp gauge	03920
Crimp gauge replacement	03919
Extra plugs and caps	03999
Rotor cap	03538

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0	50.0	PA Ultracrimp Tube	03528	25	6	27,000	79,057	29 x 108	Closure	Included	25	Plugs and Aluminum Caps
									Accessory	03538	1	Rotor cap
									Accessory	03529	1	Tube Racks
									Accessory	03920	1	Sealing Tool
50.0	43.0	Nalgene PC Oak Ridge Tube	3138-0050	50	6	27,000	79,057	29 x 108	Closure	Included	50	PP Sealing
50.0	43.0	Nalgene PPCO Oak Ridge Tube	3139-0050	50	6	27,000	79,057	29 x 107	Closure	Included	50	PP Sealing
50.0	46.0	FEP Oak Ridge Tube	3114-0050	10	6	21,500	50,129	29 x 108	Closure	Included	10	FEP Screw Top
50.0	37.0	PP Flanged Tube	03147	25	6	27,000	79,057	29 x 102	Closure	03268	25	PP Snap-on
50.0	36.0	PC Flanged Tube	03146	25	6	27,000	79,057	29 x 102	Closure	03268	25	PP Snap-on
50.0	46.0	PA Flanged Tube	03139	25	6	27,000	79,057	28 x 104	-	-	-	-
30.0	28.0	PC Oak Ridge Tube	314348	25	6	27,000	79,057	25 x 94	Closure	314347	1	Aluminum Sealing
									Adapter	75003027	2	1 place/adaptor
									Tool	314353	1	Tube Extractor Tool
30.0	28.0	PP Oak Ridge Tube	314349	25	6	27,000	79,057	25 x 94	Closure	314347	1	Aluminum Sealing
									Adapter	75003027	2	1 place/adaptor
									Tool	314353	1	Tube Extractor Tool
30.0	-	Glass Tube	-	6	6	-	-	24 x 106	Adapter	00368	1	Optional Caps
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	-	27,000	79,057	18 x 107	Closure	Included	-	PA Sealing
									Adapter	75003026	2	1 place/adaptor
16.0	16.0	Nalgene PP Oak Ridge Tube	3139-0016	50	-	27,000	79,057	18 x 107	Closure	Included	-	PA Sealing
									Adapter	75003026	2	1 place/adaptor
16.0	13.0	PP Flanged Tube	03244	50	6	27,000	79,057	18 x 100	Closure	03299	50	PP Snap-On
									Adapter	75003026	2	1 place/adaptor
15.0	-	Glass Tube	-	-	6	11,000	13,122	18 x 102	Adapter	00363	1	1 place/adaptor
14.0	11.0	PC Flanged Tube	03246	50	6	27,000	79,057	18 x 75	Closure	03269	50	PP Snap-on
									Adapter	75003025	2	1 place/adaptor
12.0	10.0	PP Flanged Tube	03116	50	6	27,000	79,057	16 x 100	Closure	03266	50	PP Nonsealing
									Adapter	75003028	2	1 place/adaptor
12.0	11.0	PC Flanged Tube	03115	50	6	27,000	79,057	16 x 100	Closure	03266	50	PP Nonsealing
									Adapter	75003028	2	1 place/adaptor
10.0	10.0	Nalgene PP Oak Ridge Tube	3139-0010	50	6	27,000	79,057	16 x 83	Closure	Included	-	PA Sealing
									Adapter	75003024	2	1 place/adaptor
10.0	9.1	PP Oak Ridge Tube ²	03929	50	6	27,000	79,057	16 x 80	Closure	03279	25	PP Sealing
									Adapter	75003024	2	1 place/adaptor
10.0	9.1	PC Oak Ridge Tube	03020	25	6	27,000	79,057	16 x 80	Closure	03279	25	PP Sealing
									Adapter	75003024	2	1 place/adaptor
10.0	10.0	Nalgene PC Oak Ridge Tube	3138-0010	50	6	27,000	79,057	16 x 83	Closure	Included	-	PA Sealing
									Adapter	75003024	2	1 place/adaptor

Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
7.0	6.0	PC Flanged Tube	03120	50	12	24,000	62,465	13 x 100	Closure	03265	50	PP Snap-On
									Adapter	00473	2	2 places/adapter
7.0	6.0	PP Flanged Tube	03121	50	12	21,000	47,825	13 x 100	Closure	03265	50	PP Snap-On
									Adapter	00473	1	2 places/adapter
4.0	3.0	PP Flanged Tube	03105	50	12	27,000	79,057	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 places/adapter
4.0	3.0	PC Flanged Tube	03104	50	12	27,000	79,057	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 places/adapter
1.5	1.5	Polyallomer Microtube	314352H01	100	18	20,300	44,826	11 x 40	Adapter	75003029	2	3 places/adapter

² Rated to 21,000 rpm when partially filled.

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ A27-8x50 Rotor

Rotor specifications

Type	Fixed angle
Material	Aluminum
Capacity	8 x 50 mL
Maximum speed	27,000 rpm
Tube dimensions Ø x L	29 x 108 mm
Angle	34°
K Factor	408
Net weight	6.8 kg/15.0 lbs
Warranty	10 years
Certifications	



	RCF	Radius
Maximum	87,207 x g	107 mm
Minimum	26,896 x g	33 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	27,000 rpm	24,000 rpm
Max. RCF	87,207 x g	68,905 x g

A27-8x50 rotor package includes:

Description	Cat. No.
A27-8x50 Rotor, including the below items:	75003008
50 mL Nalgene PPCO Oak Ridge tube with sealing cap (qty. 8)	3139-0050
Vacuum grease	76003500
Lubricant for metal parts	75003786
O-Ring replacement kit	75007009
Retaining ring pliers	65614
Instruction manual CD	50136234

Accessories:

Description	Cat. No.
Ultracrimp sealing tool and crimp gauge	03920
Ultracrimp gauge replacement	03919
Ultracrimp extra plugs and caps	03999
Rotor cap	03538

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0	50.0	PA Ultracrimp Tube	03528	25	8	27,000	87,207	29 x 108	Closure	Included	25	Plugs and Aluminum Caps
									Accessory	03538	1	Rotor cap
									Accessory	03529	2	Tube Racks
									Accessory	03920	1	Crimp Sealing Tool
50.0	37.0	PP Flanged Tube	03147	25	8	27,000	87,207	29 x 102	Closure	03268	25	PP Snap-on
50.0	36.0	PC Flanged Tube	03146	25	8	27,000	87,207	29 x 102	Closure	03268	25	PP Snap-on
50.0	46.0	PA Thin-Walled Tube	03139	25	8	21,000	52,755	28 x 104	-	-	-	-
50.0	46.0	Nalgene FEP Oak Ridge Tube	3114-0050	10	8	21,000	52,755	28 x 104	Closure	Included	10	FEP Sealing
50.0	43.0	Nalgene PC Oak Ridge Style Tube	3138-0050	50	8	27,000	87,207	29 x 107	Closure	Included	50	PP Sealing
50.0	43.0	Nalgene PPCO Oak Ridge Style Tube	3139-0050	50	8	27,000	87,207	29 x 107	Closure	Included	50	PP Sealing
30.0	28.0	Nalgene PC Oak Ridge Tube	3138-0030	50	8	27,000	87,207	26 x 102	Closure	Included	50	PP Sealing
											75003027	2
30.0	28.0	Nalgene PP Oak Ridge Tube	3139-0030	50	8	27,000	87,207	26 x 102	Closure	Included	50	PP Sealing
											75003027	2
30.0	-	Glass Tube	-	-	8	-	-	24 x 106	Adapter	00368	1	Optional Caps
16.0	16.0	Nalgene PP Oak Ridge Tube	3139-0016	50	8	27,000	87,207	18 x 107	Closure	Included	-	PA Sealing
									Adapter	75003026	2	1 place/adaptor
16.0	13.0	PP Flanged Tube	03244	50	8	27,000	87,207	18 x 100	Closure	03299	50	HDPE Sealing
									Adapter	75003026	2	1 place/adaptor
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	8	27,000	87,207	18 x 107	Closure	Included	-	PA Sealing
									Adapter	75003026	2	1 place/adaptor
15.0	15.0	Conical Tube	-	-	8	-	-	-	Adapter	75007321	2	1 place/adaptor
14.0	11.0	PC Flanged Tube	03246	50	8	27,000	87,207	18 x 75	Closure	03269	50	PP Snap-on
									Adapter	75003025	2	1 place/adaptor
12.0	10.0	PP Flanged Tube	03116	50	8	27,000	87,207	16 x 100	Closure	03266	50	PP Snap on
									Adapter	00402	2	1 place/adaptor
12.0	11.0	PC Flanged Tube	03115	50	8	27,000	87,207	16 x 100	Closure	03266	50	PP Snap on
									Adapter	00402	2	1 place/adaptor
10.0	9.1	PP Oak Ridge Tube	03929	50	8	27,000	87,207	16 x 80	Closure	03279	25	PP Snap on
									Adapter	75003024	2	1 place/adaptor
10.0	10.0	Nalgene PC Oak Ridge Tube	3138-0010	50	8	27,000	87,207	16 x 83	Closure	Included	-	PA Sealing
									Adapter	75003024	2	1 place/adaptor
10.0	10.0	Nalgene PP Oak Ridge Tube	3139-0010	50	8	27,000	87,207	16 x 83	Closure	Included	-	PA Sealing
									Adapter	75003024	2	1 place/adaptor
10.0	9.1	PC Oak Ridge Tube	03020	25	8	27,000	87,207	16 x 80	Closure	03279	25	PP Sealing
									Adapter	75003024	2	1 place/adaptor

Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
7.0	6.0	PC Flanged Tube	03120	50	8	21,000	52,755	13 x 100	Closure	03265	50	PP Snap-On
									Adapter	00473	1	1 place/adapter
7.0	6.0	PP Flanged Tube	03121	50	8	21,000	52,755	13 x 100	Closure	03265	50	PP Snap-On
									Adapter	00473	1	1 place/adapter
4.0	3.0	PP Flanged Tube	03105	50	8	27,000	87,207	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 place/adapter
4.0	3.0	PC Flanged Tubes	03104	50	16	27,000	87,207	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 place/adapter
1.5	1.5	Polyallomer Microtube	314352H01	100	24	21,000	53,302	11 x 40	Adapter	75003029	1	3 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ T29-8x50 Rotor

Rotor specifications

Type	Fixed angle
Material	Titanium
Capacity	8 x 50 mL
Maximum speed	29,000 rpm
Tube dimensions Ø x L	29 x 108 mm
Angle	34°
K Factor	354
Net weight	8.4 kg/18.5 lb
Warranty	10 years
Certifications	



	RCF	Radius
Maximum	100,605 x g	107 mm
Minimum	31,028 x g	33 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	29,000 rpm	24,000 rpm
Max. RCF	100,605 x g	68,905 x g

A29-8x50 rotor package includes:

Description	Cat. No.
T29-8x50 rotor, including the below items:	75003009
50 mL Nalgene PPCO Oak Ridge tube with sealing cap (qty. 8)	3139-0050
Vacuum grease	76003500
Lubricant for metal parts	75003786
O-Ring replacement kit	75007009
Retaining ring pliers	65614
Instruction manual CD	50136234

Accessories:

Description	Cat. No.
Ultracrimp sealing tool and crimp gauge	03920
Ultracrimp gauge replacement	03919
Extra ultracrimp plugs and caps	03999
Rotor cap	03538
Rotor stand	75003711

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0	50.0	PA Ultracrimp Tube	03528	25	8	29,000	100,605	29 x 108	Closure	Included	25	Plugs and Aluminum Caps
									Accessory	03538	1	Rotor cap
									Accessory	03529	2	Tube Racks
									Accessory	03920	1	Crimp Sealing Tool
50.0	37.0	PP Flanged Tube	03147	25	8	29,000	100,605	29 x 102	Closure	03268	25	PP Snap-on
50.0	36.0	PC Flanged Tube	03146	25	8	29,000	100,605	29 x 102	Closure	03268	25	PP Snap-on
50.0	46.0	PA Thin-Walled Tube	03139	25	8	21,000	52,755	28 x 104	-	-	-	-
50.0	46.0	Nalgene FEP Oak Ridge Tube	3114-0050	10	8	21,000	52,755	28 x 108	Closure	Included	10	FEP Sealing
50.0	43.0	Nalgene PC Oak Ridge Style Tube	3138-0050	50	8	29,000	100,605	29 x 107	Closure	Included	50	PP Sealing
50.0	43.0	Nalgene PPCO Oak Ridge Style Tube	3139-0050	50	8	29,000	100,605	29 x 107	Closure	Included	50	PP Sealing
30.0	28.0	Nalgene PC Oak Ridge Tube	3138-0030	50	8	29,000	100,605	26 x 102	Closure	Included	50	PP Sealing
									Adapter	75003027	2	1 place/adapter
30.0	28.0	Nalgene PP Oak Ridge Tube	3139-0030	50	8	29,000	100,605	26 x 102	Closure	Included	50	PP Sealing
									Adapter	75003027	2	1 place/adapter
30.0	-	Glass Tube	-	-	8	-	-	24 x 106	Adapter	00368	1	Optional Caps
16.0	16.0	Nalgene PC Oak Ridge Tube	3138-0016	50	8	29,000	100,605	-	Closure	Included	1	PA Sealing
									Adapter	75003026	2	1 place/adapter
16.0	16.0	Nalgene PP Oak Ridge Tube	3139-0016	50	8	29,000	100,605	-	Closure	Included	1	PA Sealing
									Adapter	75003026	2	1 place/adapter
16.0	13.0	PP Flanged Tube	03244	50	8	29,000	100,605	18 x 100	Closure	03299	50	HDPE Sealing
									Adapter	75003026	2	1 place/adapter
15.0	15.0	Conical Tube	-	-	8	-	-	-	Adapter	75007321	2	1 place/adapter
14.0	11.0	PC Flanged Tube	03246	50	12	29,000	100,605	18 x 75	Closure	03269	50	PP Snap-on
									Adapter	75003025	2	1 place/adapter
12.0	10.0	PP Flanged Tube	03116	50	8	29,000	100,605	16 x 100	Closure	03266	50	PP Snap on
									Adapter	75003028	2	1 place/adapter
12.0	11.0	PC Flanged Tube	03115	50	8	29,000	100,605	16 x 100	Closure	03266	50	PP Snap on
									Adapter	75003028	2	1 place/adapter
10.0	10.0	Nalgene PC Oak Ridge Tube	3138-0010	50	8	29,000	100,605	-	Closure	Included	1	PA Sealing
									Adapter	75003024	2	1 place/adapter
10.0	10.0	Nalgene PP Oak Ridge Tube	3139-0010	50	8	29,000	100,605	-	Closure	Included	1	PA Sealing
									Adapter	75003024	2	1 place/adapter
10.0	9.1	PP Oak Ridge Tube	03929	50	8	29,000	100,605	16 x 80	Closure	03279	25	PP Snap on
									Adapter	75003024	2	1 place/adapter
10.0	9.1	PC Oak Ridge Tube	03020	25	8	29,000	100,605	16 x 80	Closure	03279	25	PP Sealing
									Adapter	75003024	2	1 place/adapter


Labware ordering information continued


Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
7.0	6.0	PC Flanged Tube	03120	50	8	21,000	52,755	13 x 100	Closure	03265	50	PP Snap-On
									Adapter	00473	1	1 place/adapter
7.0	6.0	PP Flanged Tube	03121	50	8	21,000	52,755	13 x 100	Closure	03265	50	PP Snap-On
									Adapter	00473	1	1 place/adapter
4.0	3.0	PP Flanged Tube	03105	50	8	29,000	100,605	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 place/adapter
4.0	3.0	PC Flanged Tubes	03104	50	16	29,000	100,605	11 x 75	Closure	03264	50	PP Snap-On
									Adapter	75003023	2	2 place/adapter
1.5	1.5	Polyallomer Microtube	314352H01	100	24	22,600	61,492	11 x 40	Adapter	75003029	2	3 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ BIOFlex™ HC Rotor

Rotor specifications

Type	Swinging Bucket
Material	Stainless Steel with Aluminum Buckets
Capacity	4 x 1000 mL
Maximum speed	5,500 rpm
Tube dimensions Ø x L	126 x 140 mm
Angle	90°
K Factor	5,522
Net weight	5.35 kg/11.8 lbs (Rotor Cross) 10.2 kg/22.5 lbs (Cross with Buckets)
Maximum Load per Bucket	1500 g
Warranty	10 years
Certifications	

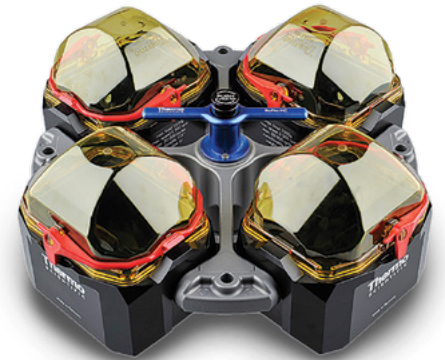

	RCF	Radius
Maximum	7,068 x g	209 mm
Minimum	3,653 x g	108 mm

	Sorvall LYNX 6000	Sorvall LYNX 4000
Max. speed	5,500 rpm	5,500 rpm
Max. RCF	7,068 x g	7,068 x g

BIOFlex™ HC rotor package includes:

Description	Cat. No.
BIOFlex™ HC rotor, including 4 buckets and including the below items:	75003000
1000 mL Polypropylene Wide-Mouth Bio-Bottles (qty. 4)	75007300
1000 mL Adapters (qty. 4)	75007301
Anti-corrosive oil	70009824
Lubricant for metal parts	75003786
Instruction manual	50136234

Accessories:

Description	Cat. No.
Thermo Scientific™ ClickSeal™ Biocontainment Lids (qty. 4)	75007309
O-Ring replacement kit for ClickSeal lids (qty. 4)	75007001
Rotor stand	75003711
Replacement BIOFlex HC buckets, set of 4	75003021
Replacement biocontainment cap for double biocontainment vessel for 50 mL conical tube (qty. 1)	50129119
Replacement O-Ring for biocontainment cap for double biocontainment vessel for 50 mL conical tube (qty. 6)	75003789
Replacement microplate carrier pad	20056846

 Biocontainment certification by Public Health England, Porton Down, UK or TÜV NORD CERT GmbH, Germany.

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
1000.0	1000.0	Polypropylene Wide-Mouth Bio-Bottle	75007300	4	4	5,500	7,068	126 x 140	Adapter	75007301	4	1 place/adapter
750.0	750.0	Polypropylene Wide-Mouth Bio-Bottle	75006443	1	4	3,600	3,028	98 x 133	Adapter	75007304	4	1 place/adapter
500.0	500.0	Fiberlite PPCO Bottle	010-1493	6	4	5,500	7,068	70 x 160	Closure	Included	6	PPGF Cap w/ PP Plug
500.0	500.0	Fiberlite PC Bottle	010-1494	6	4	5,500	7,068	70 x 160	Closure	Included	36	PPGF Cap w/ PP Plug
500.0	450.0	Nalgene PP Bottle	3141-0500	24	4	5,500	7,068	70 x 160	Adapter	75004253	4	1 place/adapter
500.0	450.0	Nalgene PC Bottle	3140-0500	24	4	5,500	7,068	70 x 160	Adapter	75004253	4	1 place/adapter
500.0	450.0	Corning™ Conical Bottle	-	-	4	5,500	7,068	98 x 133	Adapter	75007302	4	1 place/adapter
250.0	250.0	Fiberlite PPCO Bottle	010-1495	6	8	5,500	7,068	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
250.0	250.0	Fiberlite PC Bottle	010-1496	6	8	5,500	7,068	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
250.0	250.0	Nalgene PP Bottle	3141-0250	36	8	5,500	7,068	62 x 130	Adapter	75007305	4	2 places/adapter
250.0	250.0	Nalgene PC Bottle	3140-0250	36	8	5,500	7,068	62 x 130	Adapter	75007305	4	2 places/adapter
250.0	250.0	Nalgene Wide-Mouth Conical Bottle	-	-	4	5,500	7,068	26 x 145	Adapter	75005392	4	1 place/adapter
250.0	-	Corning™ Conical Bottle	-	-	4	5,500	7,068	26 x 145	Adapter	75005392	4	1 place/adapter
225.0	-	Falcon™ Conical Bottle	-	-	8	5,500	7,068	62 x 130	Adapter	75007305	4	2 places/adapter
									Adapter	BD 352090	Purchase separately	1 place/adapter
200.0	-	Nunc Conical Bottle	376813	48	8	5,500	7,068	26.5 x 139	Adapter	75007305	4	2 places/adapter
									Adapter	Nunc 377585	-	1 place/adapter
175.0	175.0	Nalgene PP Wide-Mouth Conical Bottle	3143-0175	36	8	5,500	7,068	62 x 130	Adapter	75007305	4	2 places/adapter
									Adapter	Nalgene DS3126-0175	2	1 place/adapter
175.0	175.0	Nalgene PP Wide-Mouth Conical Bottle	3144-0175	36	8	5,500	7,068	62 x 130	Adapter	75007305	4	2 places/adapter
									Adapter	Nalgene DS3126-0175	2	1 place/adapter
100.0	-	Round Bottom Open-Top Tube	-	-	30	4,500	4,732	45 x 123	Adapter	75101073	4	2 places/adapter
									BIOLink Adapter	75007304	4	-
50.0	-	Conical Tissue Culture Tube	-	-	40	5,500	7,068	29.5 x 116	Adapter	75003674	4	10 places/adapter

Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0	-	Double Biocontainment Vessels for 50 mL Conical Tube	-	-	20	5,500	7,068	29.5 x 116	Adapter	75004255	2	5 places/adapte
									Vessel	75003787	1	(10 included per 2 adapters)
		Nalgene Oak Ridge Tube	-	-	20	5,500	7,068	1 x 30	Adapter	75005802	2	1 place/adapter
		Nalgene Oak Ridge Tube	-	-	20	5,500	7,068	1 x 16	Adapter	75005803	2	1 place/adapter
		Conical Tissue Culture Tube	-	-	20	5,500	7,068	1 x 15	Adapter	75005808	2	1 place/adapter
		Nalgene Oak Ridge Tube	-	-	20	5,500	7,068	1 x 16	Adapter	75005803	2	1 place/adapter
		Blood Collection Tube	-	-	20	5,500	7,068	1 x 10	Adapter	75005804	2	1 place/adapter
		Blood Collection Tube	-	-	20	5,500	7,068	1 x 7	Adapter	75005805	2	1 place/adapter
		Blood Collection Tube (11)	-	-	40	5,500	7,068	2 x 3.5	Adapter	75005806	2	1 place/adapter
		Microtube	-	-	40	5,500	7,068	2 x 1.5/2	Adapter	75005807	2	1 place/adapter
50.0	42.0	Nalgene PP Oak Ridge Tube	3139-0050	100	48	5,500	7,068	29.5 x 120	Adapter	75004252	4	1 place/adapter
50.0	42.0	Nalgene PP Oak Ridge Tube	3138-0050	100	48	5,500	7,068	29.5 x 120	Adapter	75004252	4	12 places/adapter
50.0	50.0	Conical or Skirted Tube	-	-	20	5,500	7,068	29.5 x 120	Adapter	75003824	4	12 places/adapter
									BIOLink Adapter	75007304	4	-
25.0	25.0	Universal Conical or Skirted Tube	-	-	28	4,500	4,732	29.5 x 120	Adapter	75003716	4	7 places/adapter
									BIOLink Adapter	75007304	4	-
25.0	25.0	Universal Conical or Skirted Tube	-	-	28	4,500	4,732	29.5 x 120	Adapter	75003716	4	7 places/adapter
									BIOLink Adapter	75007304	4	-
16.0	16.0	Nalgene PP Round Bottom Tube	3139-0016	50	48	4,500	4,732	18 x 134	Adapter	75003718	4	12 places/adapter
									BIOLink Adapter	75007304	4	-
16.0	16.0	Nalgene PC Round Bottom Tube	3138-0016	50	48	4,500	4,732	18 x 134	Adapter	75003718	4	12 places/adapter
									BIOLink Adapter	75007304	4	-
15.0	-	Conical Tissue Culture Tube	-	-	96	5,500	7,068	17.5 x 121	Adapter	75007306	4	24 places/adapter

Labware ordering information continued


Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./ pkg.	Tubes/ rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
15.0	-	Blood Collection Tube	-	-	32	4,500	4,732	17 x 125	Adapter	75003719	4	8 places/adapter
									BIOLink Adapter	75007304	4	-
10	-	Blood Collection Tube or Corex™/ Kimble™ Tube	-	-	148	4,500	4,732	17 x 110	Adapter	75003672	4	37 places/adapter
3	-	RIA or Round Bottom Tube (without cap)	-	-	148	4,500	4,732	13 x 116	Adapter	75003724	4	37 places/adapter
									BIOLink Adapter	75007304	4	-
4.5/6 ml		Blood Collection Tube (Greiner™)	-	-	164	5,500	7,068	14 x 110	Adapter	75003709	4	41 places/adapter
5/7 ml		Blood Collection Tube (BD)	-	-	196	4,500	4,732	14 x 110	Adapter	75003671	4	49 places/adapter
1.5 - 2.0		Conical Microtube	-	-	192	5,500	7,068		Adapter	75003733	4	48 places/adapter
									BIOLink Adapter	75007304	4	-
Microplates		Standard Microplate	-	-	24	5,500	7,068	86 x 128	Adapter	75007303		6 places/adapter
Microplates		Deep well Microplate	-	-	8	5,500	7,068	86 x 128	Adapter	75007303		2 places/adapter
Flasks		T-75 Nunc Easy Flask	-	-	4	2,925	2,000	11 x 50	Adapter	75008383	4	1 place/adapter
									BIOLink Adapter	75007304	4	-
Flasks		T-25 Nunc Easy Flask	-	-	8	2,925	2,000		Adapter	75008384	4	1 place/adapter
									BIOLink Adapter	75007304	4	-
Blood Bags	-	Small Blood Bag / Cell Culture Bags	-	-	-	3,600	3,028		Adapter	75003829	4	2 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ BIOFlex™ HS Rotor



Rotor specifications

Type	Swinging Bucket
Material	Stainless Steel with Aluminum Buckets
Capacity	4 x 1000 mL
Maximum speed	7,000 rpm
Tube dimensions Ø x L	80 x 125 mm
Angle	90°
K Factor	4,889
Net weight	4.36 kg/9.6 lbs (Rotor Cross) 7.64 kg/16.8 lbs (Cross with Buckets)
Maximum Load per Bucket	600 g
Warranty	10 years
Certifications	


	RCF	Radius
Maximum	10,025 x g	183 mm
Minimum	3,890 x g	71 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	7,000 rpm	7,000 rpm
Max. RCF	10,025 x g	10,025 x g

BIOFlex™ HS rotor package includes:

Description	Cat. No.
BIOFlex™ HS rotor, including 4 Buckets and including the below items:	75003002
400 mL Polypropylene Wide-Mouth Bio-Bottles (qty. 4)	75007585
Anti-corrosive oil	70009824
Lubricant for metal parts	75003786
Instruction manual	50136234

Accessories:

Description	Cat. No.
ClickSeal biocontainment lids (qty. 4) 	75003656
O-Ring replacement kit for ClickSeal lids (qty. 4)	75003657
Rotor stand	75003711
Replacement BIOFlex HS buckets, set of 4	75003040
Replacement O-Ring for biocontainment Cap for double biocontainment vessel for 50 mL conical tube (qty. 6)	75003789

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
400.0	400.0	Polypropylene Bio-Bottle	75007585	12	4	7,000	10,025	80 x 125	-	-	-	-
250.0	250.0	Fiberlite PPCO Bottle	010-1495	6	4	7,000	10,025	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
250.0	250.0	Fiberlite PC Bottle	010-1496	6	4	7,000	10,025	61 x 124	Closure	Included	6	PPGF Cap w/ PP Plug
250.0	250.0	Nunc Wide-Mouth Conical Bottle	376814	40	4	7,000	10,025	60x144	Adapter	75004258	4	1 place/adapter
250.0	250.0	Nalgene PP Bottle	3141-0250	36	4	7,000	10,025	62 x 135	Adapter	75004257	4	1 place/adapter
250.0	250.0	Nalgene PC Bottle	3140-0250	36	4	7,000	10,025	62 x 135	Adapter	75004257	4	1 place/adapter
225.0	-	Falcon™ Conical Bottle	-	-	4	4,800	4,700	62 x 130	Adapter	75004257	4	1 place/adapter
									Adapter	BD 352090	Purchase Separately	1 place/adapter
200.0	200.0	Nunc PP Conical Bottle	376813	48	4	7,000	10,025	62 x 125	Adapter	75004258	4	1 place/adapter
175.0	175.0	Nalgene PP Wide-Mouth Conical Bottle	3143-0175	36	4	7,000	10,025	62 x 125	Adapter	75004258	4	1 place/adapter
175.0	175.0	Nalgene PC Wide-Mouth Conical Bottle	3144-0175	36	4	7,000	10,025	62 x 125	Adapter	75004258	4	1 place/adapter
100.0	-	Round Bottom Open-Top Tube	-	-	4	4,800	4,700	45 x 117	Adapter	75003708	4	1 place/adapter
50.0	-	DIN Round Bottom Tube	-	-	12	4,800	4,700	34.5 x 105	Adapter	75003707	4	3 places/adapter
50.0	42.0	Nalgene PP Oak Ridge Tube	3139-0050	100	16	7,000	10,025	28.5 x 114	Adapter	75003799	4	4 places/adapter
50.0	42.0	Nalgene PC Oak Ridge Tube	3138-0050	100	16	7,000	10,025	28.5 x 114	Adapter	75003799	4	4 places/adapter
50.0	-	Conical Tissue Culture Tube	-	-	16	4,800	4,700	29.5 x 116	Adapter	75003683	4	4 places/adapter
50.0	-	Conical Tissue Culture Tube	-	-	12	7,000	10,025	29.5 x 116	Adapter	75005393	4	3 places/adapter
50.0	-	Double Biocontainment Vessels for 50 mL Conical Tube	-	-	12	7,000	10,025	29.5 x 116	Adapter	75004259	2	2 places/adapter
									Vessel	75003787	1	(2 included per 2 adapters)
		Nalgene Oak Ridge Tube	-	-	12	7,000	10,025	1 x 30	Adapter	75005802	2	1 place/adapter
		Nalgene Oak Ridge Tube	-	-	12	7,000	10,025	1 x 16	Adapter	75005803	2	1 place/adapter
		Conical Tissue Culture Tube	-	-	12	7,000	10,025	1 x 15	Adapter	75005808	2	1 place/adapter
		Blood Collection Tube	-	-	12	7,000	10,025	1 x 10	Adapter	75005804	2	1 place/adapter
		Blood Collection Tube	-	-	12	7,000	10,025	1 x 7	Adapter	75005805	2	1 place/adapter
		Blood Collection Tube (11 mm)	-	-	24	7,000	10,025	2 x 3.5	Adapter	75005806	2	2 places/adapter
		Microtube	-	-	24	7,000	10,025	2 x 1.5/2	Adapter	75005807	2	2 places/adapter


Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
30.0	-	DIN Round/Flat Bottom Tube	-	-	20	7,000	10,025	25.5 x 108	Adapter	75003703	4	5 places/adapter
20.0	-	Round Bottom Tube	-	-	12	4,800	4,700	25 x 110	Adapter	75003706	4	3 places/adapter
16.0	16.0	Nalgene PP Round Bottom Tube	3139-0016	50	28	7,000	10,025	18 x 112	Adapter	75003798	4	7 places/adapter
16.0	16.0	Nalgene PC Round Bottom Tube	3138-0016	50	28	7,000	10,025	18 x 112	Adapter	75003798	4	7 places/adapter
15.0	-	Conical Tissue Culture Tube	-	-	36	7,000	10,025	17 x 121	Adapter	75005394	4	9 places/adapter
15.0	-	Blood Collection Tube (17 x 125 mm)	-	-	16	4,800	4,700	15.5 x 131	Adapter	75003794	4	4 places/adapter
15.0	-	Round Bottom Tube (Sarstedt™)	-	-	40	4,800	4,700	17 x 105	Adapter	75003704	4	10 places/adapter
10/15 ml	-	Corex™/Kimble™ Tube or 10 mL Blood Collection (BD Vacutainer™/Vacuette™)	-	-	56	7,000	10,025	17 x 113	Adapter	75003681	4	14 places/adapter
5/7 ml	-	Blood Collection Tube (Vacutainer)	-	-	76	4,800	4,700	13 x 110	Adapter	75003680	4	19 places/adapter
4.5/6 ml	-	Blood Collection Tube (Greiner™)	-	-	64	7,000	10,025	13 x 110	Adapter	75003825	4	16 places/adapter
3/5 ml	-	RIA or Round Bottom Tube (without cap)	-	-	76	7,000	10,025	11 x 110	Adapter	75003793	4	19 places/adapter
1.5-2ml	-	Conical/Round Microtube	-	-	136	7,000	10,025	11 x 45	Adapter	75003700	4	34 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ TH13-6x50 Rotor

Rotor specifications

Type	Swinging Bucket
Material	Aluminum with Titanium Buckets
Capacity	6 x 50 mL
Maximum speed	13,100 rpm
Tube dimensions Ø x L	29 x 104 mm
Angle	90°
K Factor	1,503
Net weight	7.2 kg/15.8 lbs
Warranty	10 years
Certifications	



	RCF	Radius
Maximum	30,314 x g	158 mm
Minimum	10,936 x g	57 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	13,100 rpm	13,100 rpm
Max. RCF	30,314 x g	30,314 x g

TH13-6x50 rotor package includes:

Description	Cat. No.
TH13-6x50 rotor with buckets and biocontainment covers and including the below items:	75003010
Rotor cleaning kit	12259
Instruction manual CD	501362 34
Vacuum grease	76003500
Lubricant for metal parts	75003786

Accessories:

Description	Cat. No.
Replacement biocontainment covers (each)	50129119
Bucket O-Rings replacement, set of 6	75007002

Labware ordering information

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
50.0		PA Thin-Walled Tube	03139	25	6	13,100	30,314	28 x 104	-	-	-	-
50.0	50.0	Nunc PP Conical Tube	339650	25	6	-	-	-	Adapter	75004264	1	1 place/adapter
50.0	50.0	Falcon™ PP Conical Tube	-	-	6	13,100	30,314	-	Adapter	75004264	1	1 place/adapter
50.0	50.0	Corning™ Conical Tube	-	-	6	13,100	30,314	-	Adapter	75004264	1	1 place/adapter
50.0	50.0	Sarstedt™ Conical Tube	-	-	6	13,100	30,314	-	Adapter	75004264	1	1 place/adapter
50.0	50.0	Sterilin™ PP Conical Tube	-	-	6	13,100	30,314	-	Adapter	75004264	1	1 place/adapter
50.0	50.0	Conical Filter Tube (e.g. Amicon™)	-	-	6	-	-	-	Adapter	75004264	1	1 place/adapter
50.0	50.0	Greiner™ Conical Tube	-	-	6	13,100	30,314	-	Adapter	75004264	1	1 place/adapter
50.0	48.0	PC Flanged Tube	03146	25	6	13,100	30,314	29 x 102	Closure	03268	25	PP Snap-on
50.0	46.0	PP Flanged Tube	03147	25	6	13,100	30,314	29 x 102	Closure	03268	25	PP Snap-on
50.0	43.0	Nalgene PC Oak Ridge Tube	3118-0050	100	6	13,100	30,314	29 x 107	Closure	Included	100	PP Screw Top
50.0	43.0	Nalgene POCO Oak Ridge Tube	3119-0050	50	6	13,100	30,314	29 x 107	Closure	Included	100	PP Screw Top
50.0	43.0	Nalgene PC Oak Ridge Tube	3138-0050	50	6	13,100	30,314	29 x 107	Closure	Included	50	PP Sealing
50.0	43.0	Nalgene POCO Oak Ridge Tube	3139-0050	50	6	13,100	30,314	29 x 107	Closure	Included	50	PP Sealing
16.0	16.0	PP Flanged Tube	03244	50	6	13,100	30,314	18 x 100	Closure	03299	50	PP Snap-on
									Adapter	75003026	2	1 place/adapter
15.0	15.0	Nunc PP Conical Tube	339650	50	6	-	-	-	Adapter	75007321	2	1 place/adapter
15.0	15.0	Greiner™ Conical Tube	-	-	6	13,100	30,314	-	Adapter	75007321	2	1 place/adapter
15.0	15.0	Falcon™ PP Conical Tube	-	-	6	13,100	30,314	-	Adapter	75007321	2	1 place/adapter
15.0	15.0	Corning™ Conical Tube	-	-	6	13,100	30,314	-	Adapter	75007321	2	1 place/adapter
15.0	15.0	Sarstedt PP Conical Tube	-	-	6	13,100	30,314	-	Adapter	75007321	2	1 place/adapter
15.0	15.0	Sterilin™ PP Conical Tube	-	-	6	13,100	30,314	-	Adapter	75007321	2	1 place/adapter
15.0	15.0	Conical Filter Tube (e.g. Amicon™)	-	-	6	-	-	-	Adapter	75007321	2	1 place/adapter
14.0	11.0	PC Flanged Tube	03246	50	12	13,100	30,314	18 x 75	Closure	03269	50	PP Snap-on
									Adapter	75003025	2	1 place/adapter
12.0	12.0	PP Flanged Tube	03116	50	6	13,100	30,314	16 x 100	Closure	03266	50	PP Snap-on
									Adapter	75003028	2	1 place/adapter
12.0	12.0	PC Flanged Tube	03115	50	6	13,100	30,314	16 x 100	Closure	75003028	50	PP Snap-on
									Adapter	00367	2	1 place/adapter

Labware ordering information continued

Tube vol. (mL)	Vol. fill (mL)	Description	Cat. No.	Qty./pkg.	Tubes/rotor	Max. speed (rpm)	Max. RCF ² (x g)	Dims ØxL (mm)	Required tube closures, adapters, tools & accessories			
									Type	Cat. No.	Qty./pkg.	Description
10.0	-	Pyrex™ Graduated Conical Tube	-	-	6	5,000	4,416	18 x 100	Adapter	00367	1	1 place/adapter
10.0	9.7	PC Oak Ridge Tube	03020	50	6	13,100	30,314	16 x 83	Closure	03279	25	PP Sealing
									Closure	03924	25	PP Screw Top
									Adapter	75003024	2	1 place/adapter
10.0	9.1	PP Oak Ridge Tube	03929	50	6	13,100	30,314	16 x 83	Closure	03279	25	PP Sealing
									Closure	03924	25	PP Screw Top
									Adapter	75003024	2	1 place/adapter
7.0	7.0	PC Flanged Tube	03120	50	6	13,100	30,314	13 x 100	Closure	03265	50	PP Snap-on
									Adapter	00473	1	1 place/adapter
7.0	7.0	PP Flanged Tube	03121	50	6	13,100	30,314	13 x 100	Closure	03265	50	PP Snap-on
									Adapter	00473	1	1 place/adapter
4.0	4.0	PP Flanged Tube	03105	50	12	13,100	30,314	11 x 75	Closure	03264	50	PP Snap-on
									Adapter	00473	1	1 place/adapter
4.0	4.0	PC Flanged Tube	03104	50	12	13,100	30,314	11 x 75	Closure	03264	50	PP Snap-on
									Adapter	75003023	2	2 places/adapter
1.5	1.5	Polyallomer Microtube	314352H01	100	18	10,100	18,112	11 x 40	Adapter	75003029	2	3 places/adapter
1.0	1.0	Cellulose (Acetate Butyrate) Tube	03103	50	24	13,100	30,314	7 x 50	Adapter	00408	1	4 places/adapter

Refer to manufacturer's recommendations for product performance information.

Thermo Scientific™ TCF-20 Continuous Flow and Zonal Rotors



Rotor specifications	
Type	Continuous Flow / Zonal
Material	Titanium
Capacity	1350 mL
Maximum speed	20,000 rpm
K Factor	620
Net weight	7.8 kg/17.2 lbs
Warranty	5 years

	RCF	Radius
Maximum	42,931 x g	96.0 mm
Minimum	16,099 x g	36.0 mm

	Sorvall LYNX 6000 Centrifuge	Sorvall LYNX 4000 Centrifuge
Max. speed	20,000 rpm	N/A
Max. RCF	42,931 x g	N/A

TCF-20 continuous flow rotor package includes:

Description	Cat. No.
TCF-20 continuous flow rotor package including the below items:	75003012 #
Titanium rotor body	75000604
Rotor lid	75004012
Septa	49609
Tube holder assembly	75000601
Tool kit	75004052
Continuous flow assembly	75000605
Instruction manual	50139705

TCF-20 zonal flow rotor package includes:

Description	Cat. No.
TCF-20 zonal flow rotor package	75003013 #
Titanium rotor body	75000604
Rotor lid	75004012
Septa	49609
Tube holder assembly	75000601
Tool kit	75004052
Zonal mode assembly	75004051
Instruction manual	50139705

Continuous flow kit 75007637 is needed to order separately for rotor operation.

Accessories:

Description	Cat. No.
Peristaltic pump	20210448
Latex tubing	13117
O-Ring set for continuous flow assembly	75004006
O-Ring set for zonal mode assembly	75004007
Screws and small parts	75004008
Bearing kit	75004050
Tube holder assembly	75000601
Titanium rotor body	75000604
Rotor lid	75004012
Septa	49609
Tool kit	75004052
Continuous flow assembly	75000605
Zonal mode assembly	75004051
Rotor unlocking tool	75000602

 Learn more at thermofisher.com/lynx

For Laboratory Use. It is the customer's responsibility to ensure that the performance of the product is suitable for customers' specific uses or applications. © 2023 Thermo Fisher Scientific Inc. All rights reserved. Falcon™ and Corning™ are registered trademarks of Corning™. BD Vacutainer is a registered trademark of Becton, Dickinson and Company. All other trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details. **EXT6117 1123**